

560 Series DISS Medical Gas Outlets

Elements

- Safety keyed to help prevent cross-connections
- Conforms to NFPA 99 and CGA/DISS standards
- Modular design
- Cleaned for oxygen service
- UL listed



Description

The Chemetron® 560 Series DISS outlet stations provide convenient access to medical gases in patient care locations. They are fully compatible with DISS style diameter indexed threaded adapters.

The outlets are field-gangable in any order desired by the facility. Outlets are spaced 5" (12.7 cm) center to center, and may be ganged with Chemetron modular electrical assemblies in a console arrangement. Outlet finish plates and rough-in plates are pin indexed using steel pins for positive cross-connection prevention. Primary valves are made of brass, chrome-plated, and contain a self-sealing poppet. Primary valves are field adjustable to compensate for variations in plaster thickness of 1/2" to 3/4". Outlets may be special ordered with extended valve bodies that allow for an additional 3/8" or 3/4" adjustment.

The type K copper inlet tube has a 1/2" OD connection, and can rotate 360 degrees to accommodate any field piping arrangement.



560 Series DISS Medical Gas Outlet

The primary gas check can be serviced without shutting off the gas supply.


The outlet stations are fully compliant with NFPA 99 1999, CGA V-5 and CSA Z-305, and are cleaned for oxygen service in accordance with CGA G-4.1. The outlet is UL listed. Rough-in plates for the outlets are zinc plated steel for corrosion resistance and fire safety. Finish plates are die-cast zinc, chrome-plated with brushed stainless dress plate for appearance. Gas identification is provided by a polycarbonate plastic insert, color-coded by gas service.

Vacuum outlets may be used with or without slide brackets for vacuum collection bottles of up to 1/2-gallon capacity. (Where slide brackets are used, order one per vacuum outlet.)

Ordering Information

Gas Service	Catalog Number U.S. Color Codes	Quantity	Catalog Number ISO Color Codes	Quantity
Oxygen	64-04-5601		64-04-5651	
Vacuum	64-04-5602		64-04-5652	
Air	64-04-5603		64-04-5653	
Nitrous Oxide	64-04-5604		64-04-5604	
Evacuation Vacuum	64-04-5606		64-04-5606	
Carbon Dioxide	64-04-5607		64-04-5607	
Oxygen/CO2 (CO2 less than 7%)	64-04-5608		64-04-5658	
Slide Bracket	64-06-0001		64-06-0001	
Nitrogen	64-04-5605		N/A	

Dimensional Data

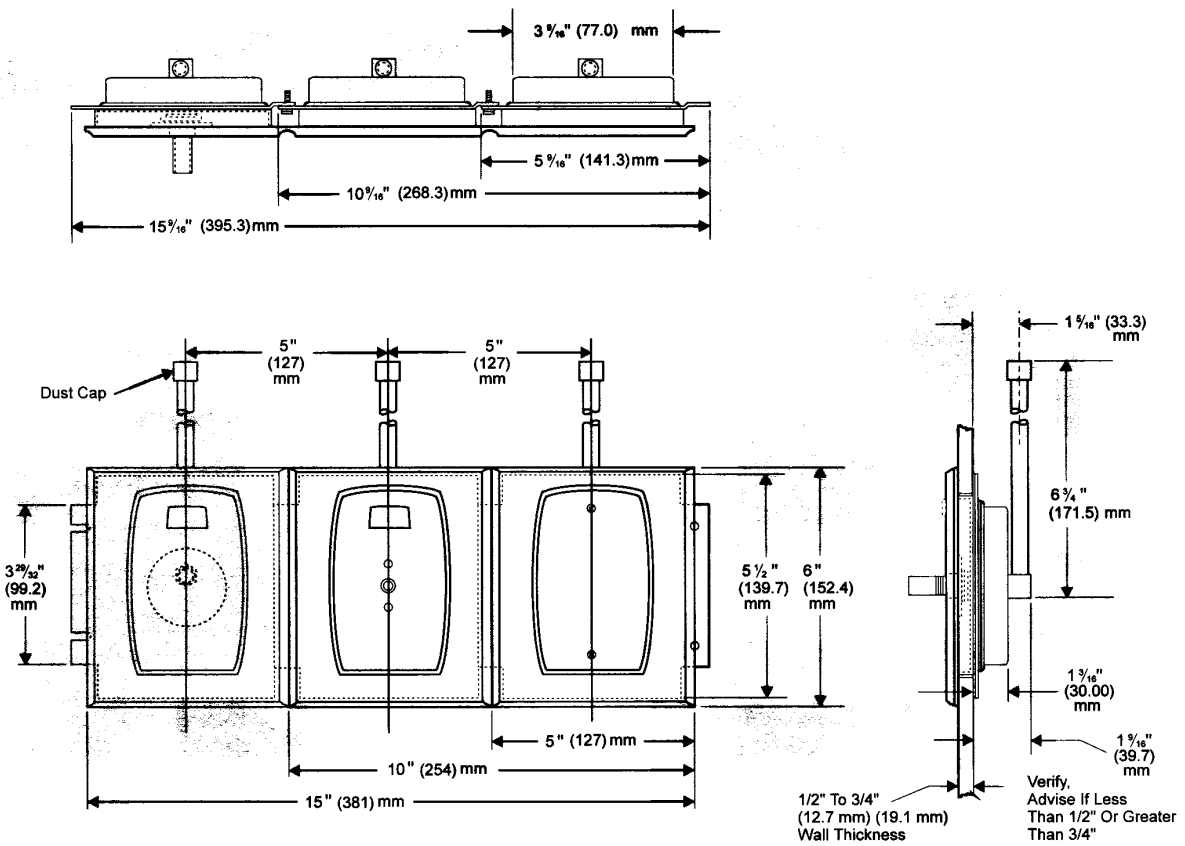

Refer to installation instructions prior to assembly.



560 Series DISS Medical Gas Outlets

All specifications are nominal and subject to change without notice.

Warranty:
See Allied Statement of Warranties for details.

© 2005 Allied Healthcare Products, Inc.® Printed in USA

Customer Service:

USA: (800) 444-3960

FAX: (314) 771-7537

Canada/Mexico: (800) 446-0552

FAX: (800) 246-6201

International: (314) 268-1683

FAX: (314) 771-5183



1720 Sublette Avenue
St. Louis, MO 63110 USA
(314)771-2400

www.alliedhpi.com

ISO 9001 CERTIFIED

Product Specification